


Pathfinder[®] is well suited for:

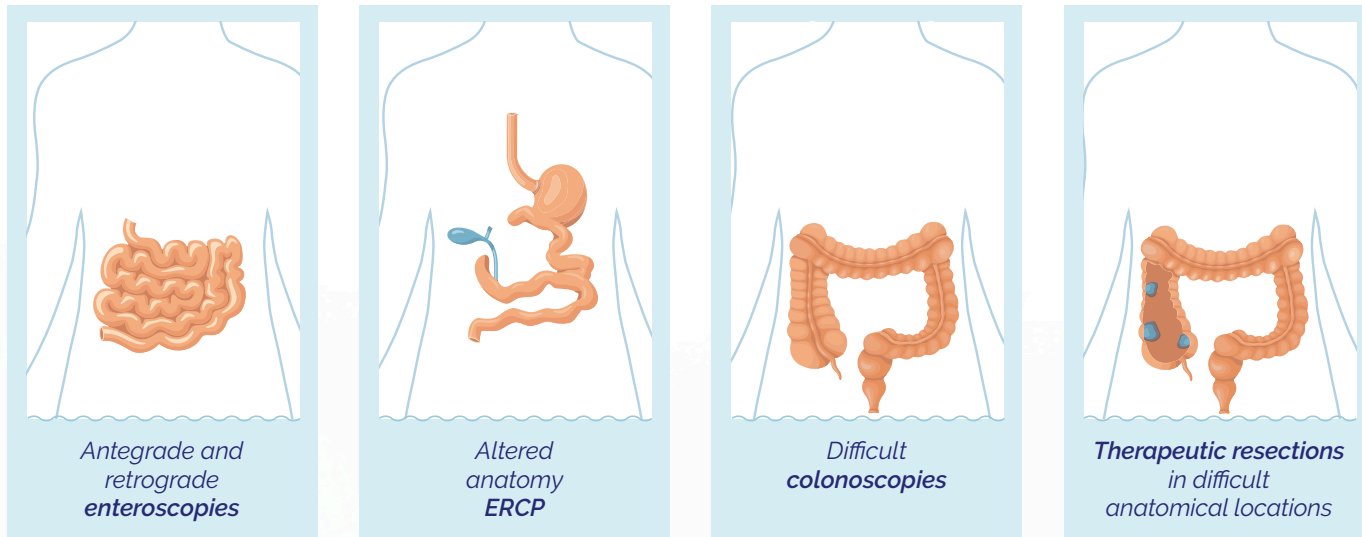

Single use device.
Do not reuse.
Do not attempt to resterilize.


PATHFINDER Ordering Information				
Catalog No	Scope Compatibility	Inner Diameter	Outer Diameter	Total Length
GI 085140-2	11.5 mm-11.8 mm	14 mm	16.5 mm	85 cm
GI 110140-2				110 cm
GI 085160-2	12.8 mm	16 mm	18.5 mm	85 cm
GI 085160-2				110 cm

PATHFINDER ENDOSCOPE CAP Ordering Information				
Catalog No	Scope Compatibility	Max Insertion Diameter (OD)	Depth	Total Length
GI EC 0412-1	11.7 mm tip diameter	12.2 mm	4 mm	8 mm

References

1: Test data on file at Neptune Medical.

FINALLY...
AN EASY WAY TO
TAME LOOPING

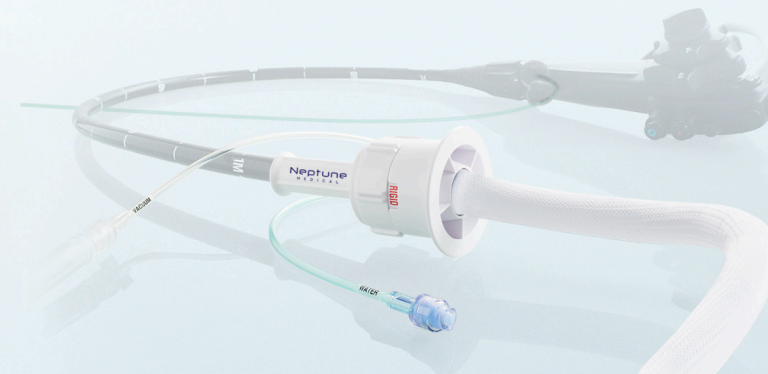

Precisely Control
Scope Tip

Reduce Looping
Recurrence

Easily Swap or
Reinsert Scope



Multiple U.S. and foreign patents pending.
APM006/E/DCO-0177
© 2021 Neptune Medical

Neptune Medical
1828 El Camino Real, Suite 508
Burlingame, CA 94010-3118

info@neptunemedical.com
888-202-8401
www.GIPathfinder.com

Caution: This device is cleared for sale only in the U.S. Federal law restricts sale by or on the order of a physician. Consult the information for use (IFU) insert provided with the product for indications, contraindications, warnings, precautions and potential adverse events.

Pathfinder® is a novel rigidizing overtube designed to manage loop formation, help you easily navigate to your target, and control your endoscope tip for diagnostic and therapeutic GI procedures.

CLINICAL BENEFITS

Mitigate looping recurrence
Goes from flexible to rigid and back again—quickly, easily and repeatedly—to prevent looping from complicating your procedures.

Precisely control endoscope tip
Achieve one-to-one motion in almost any anatomy to allow your scope tip to move and torque exactly as you intend, making challenging maneuvers faster & easier.

Easily swap or reinsert scope in-place
Maintains an open lumen without ovalization and kinks, allowing scope swap and enabling easy removal of large specimens or multiple polyps.

Make procedures easier for you and your team
Virtually eliminates the need for abdominal pressure and position changes to manage looping, reducing frustration and exertion.

Extremely Flexible
Matches curves of the endoscope

Atraumatic Tip
Minimal tip gap enhances safety

Durable Hydrophilic Coating
Sustained benefits for long procedures

Proximal Seal
Prevents back flow during flush and maintains insufflation

Endoscope Cap
Proprietary thin-walled endoscope cap allows easy swapping of PCF colonoscopes

WATCH
PATHFINDER
IN ACTION



Scan code with your phone camera to view video



ABOUT DYNAMIC RIGIDIZATION™

Vacuum control activates thin, microcomposite layers in overtube walls to go from flexible to rigid, and back again, on-demand.



Intuitively toggle between Flex & Rigid states



No foreshortening or shape change

15x

Up to 15x change in stiffness¹



Length or tortuosity do not degrade stiffness



No expensive capital equipment



Single use device